

chamber accessories

Slice Anchors for Series 20 and 40 Chambers

Provides ability to use tissue sections in most RC Series chambers



Slice anchors are available for most open bath chambers in the 20 and 40 Series family except for field stimulation chambers. Anchors are constructed for an easy push-in fit into the chamber's bath area. This allows the user to control the anchor pressure applied to slice. Most anchors are made from type 316 stainless steel and Lycra® and finished with a plastic coating. Plastic anchors are available for the RC-27L Chamber.

Order #	Model	Product
W4 64-0246	SHD-22L/15	For RC-22 Chamber, 1.5 mm
W4 64-0247	SHD-22CL/15	For RC-22C Chamber, 1.5 mm
W4 64-0248	SHD-22L/10	For RC-22 Chamber, 1.0 mm
W4 64-0249	SHD-22CL/10	For RC-22C Chamber, 1.0 mm
W4 64-0250	SHD-26H/10	For RC-26 Chamber, 1.0 mm
W4 64-0251	SHD-26H/15	For RC-26 Chamber, 1.5 mm
W4 64-0252	SHD-26H/2	For RC-26 Chamber, 2.0 mm
W4 64-0253	SHD-26GH/10	For RC-26G Chamber, 1.0 mm
W4 64-0254	SHD-26GH/15	For RC-26G Chamber, 1.5 mm
W4 64-0255	SHD-26GH/2	For RC-26G Chamber, 2.0 mm
W4 64-0256	SHD-27LH/10	For RC-27L & RC-29, 1.0 mm
W4 64-0257	SHD-27LH/15	For RC-27L & RC-29, 1.5 mm
W4 64-0258	SHD-27LH/2	For RC-27L & RC-29, 2.0 mm
W4 64-0194	SHD-27H/10	For RC-27 Chamber, 1.0 mm
W4 64-0193	SHD-27H/15	For RC-27 Chamber, 1.5 mm
W4 64-0192	SHD-27H/2	For RC-27 Chamber, 2.0 mm
W4 64-0197	SHD-27N/10	For RC-27N Chamber, 1.0 mm

Order #	Model	Product
W4 64-0196	SHD-27N/15	For RC-27N Chamber, 1.5 mm
W4 64-0195	SHD-27N/20	For RC-27N Chamber, 2.0 mm
W4 64-1411	3SHD-22F/15	Flat for RC-22 Chamber, 1.5 mm
W4 64-1412	SHD-22F/10	Flat for RC-22 Chamber, 1.0 mm
W4 64-1413	SHD-22CF/15	Flat for RC-22C Chamber, 1.5 mm
W4 64-1414	SHD-22CF/10	Flat for RC-22C Chamber, 1.0 mm
W4 64-1415	SHD-40/2	Flat for RC-40 Chamber, 2.0 mm
W4 64-1416	SHD-40/15	Flat for RC-40 Chamber, 1.5 mm
W4 64-1417	SHD-40/10	Flat for RC-40 Chamber, 1.0 mm
W4 64-1418	SHD-41/15	Flat for RC-41 Chamber, 1.5 mm
W4 64-1419	SHD-41/10	Flat for RC-41 Chamber, 1.0 mm
W4 64-1420	SHD-42/15	Flat for RC-42 Chamber, 1.5 mm
W4 64-1421	SHD-42/10	Flat for RC-42 Chamber, 1.0 mm
W4 64-0259	SHD-27LP/10	Plastic for RC-27L, RC-27LD & RC-29, 1.0 mm
W4 64-0260	SHD-27LP/15	Plastic for RC-27L, RC-27LD & RC-29, 1.5 mm
W4 64-0261	SHD-27LP/2	Plastic for RC-27L, RC-27LD & RC-29, 2.0 mm

Slice Anchor Kits and Supports for Series 20 Chambers

Provides ability to use tissue sections in most RC Series chambers

Slice Anchor Kits

Slice Anchor Kits are available for all Series 20 tissue slice recording chambers. Kits for RC-26, RC-26G, RC-27, and RC-27L, (large bath chambers) contain three anchors with 0.1 mm thick Lycra® threads at 1.0, 1.5, and 2.0 mm spacing. Kits are available with plastic anchors or stainless steel.

Kits for the RC-22 and RC-22C (smaller slice chambers) contain two stainless steel anchors with Lycra® threads at 1.0 mm and 1.5 mm spacing. The stainless steel material allows the anchors to be adjusted for a closer or looser fit in chamber.

All kits come with nylon mesh in two sizes (112 micron and 160 micron) that may be cut and used in conjunction with the slice anchor.




Order #	Model	Product
W4 64-0263	SHD-22KIT	Slice Anchor Kit for RC-22
W4 64-0264	SHD-22CKIT	Slice Anchor Kit for RC-22C
W4 64-0265	SHD-26KIT	Slice Anchor Kit for RC-26
W4 64-0266	SHD-26GKIT	Slice Anchor Kit for RC-26G
W4 64-0267	SHD-27KIT	Slice Anchor Kit for RC-27
W4 64-0268	SHD-27LKIT	Slice Anchor Kit for RC-27L and RC-29
W4 64-0190	SHD-27LPKIT	Plastic Slice Anchor Kit for RC-27L and RC-29 Chambers
W4 64-0191	SHD-27NKIT	Slice Anchor Kit for RC-27N

Slice Support with Nylon Grid

The SS-3G is a version of our standard SS-3 providing additional support by replacing the Lycra strands with a nylon mesh. The support is machined from 0.5 mm polycarbonate and the mesh is 112 micron x 0.1 micron thick.

The SS-4, SS-4-500V and SS-4-500H with Lycra thread provide additional slice support by using strands on a narrower 500 micron spacing. The support is machined from 0.5 mm polycarbonate and the Lycra strands are available in either horizontal or vertical configurations.

Order #	Model	Product
W4 64-1533	SS-3G	Slice Support for RC-27L Chamber with grid
W4 64-1550	SS-4-500V	Vertical Slice Support for RC-27D and RC-27LD Drip Chambers
W4 64-1551	SS-4-500H	Horizontal Slice Support for RC-27D and RC-27LD Drip Chambers

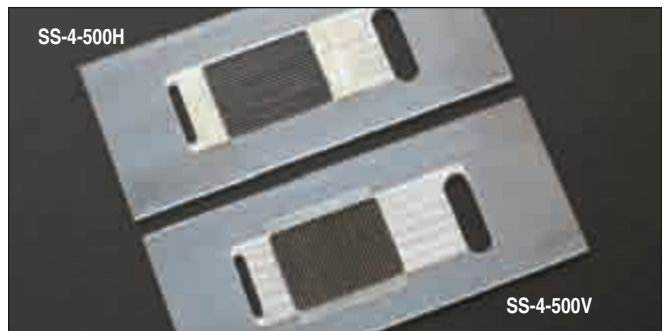

Nylon Mesh Kit

Ideal for supporting and holding slices in recording and imaging chambers, this nylon mesh can be used in conjunction with Warner slice anchors or other holders. Each kit comes with nylon mesh in two sizes 0.3 mm and 1.13 mm mesh opening and are 0.2 and 0.12 mm thick respectively. Five 7 x 7 cm squares of each size are supplied.

Order #	Model	Product
W4 64-0198	NYL/MESH	Nylon Mesh Kit – 10 pieces

